Hello Bio, Inc. 304 Wall St., Princeton, NJ 08540 USA

T. 609-683-7500 F. 609-228-4994

customercare-usa@m2stage.hellobio.com



# **DATASHEET**

Picro Sirius Red Solution

# **Product overview**

Name Picro Sirius Red Solution

Cat No HB9475

Alternative names PSR Stain Solution, Picrosirius Red Solution

Biological description Overview

The picro sirius red staining technique is frequently used to histologically stain collagen I and III fibers. It can be used to identify fibrillar collagen networks in tissue sections.

When viewed under linear polarized light, collagens appear as green, red, orange or yellow fibers.


Sample orientation under linear polarized light is thought to affect the hue and signal strength of the technique. Circular polarized light may be used to overcome this limitation but requires specialized equipment. Alternatively, fluorescent imaging of picro sirius red stained samples may be used which yields a strong red fluorescent signal that is sensitive and specific for collagen and is unaffected by


sample orientation. Dyes & stains

**Description** Solution for histological staining of Collagen I and III fibers

# **Images**

**Biological action** 

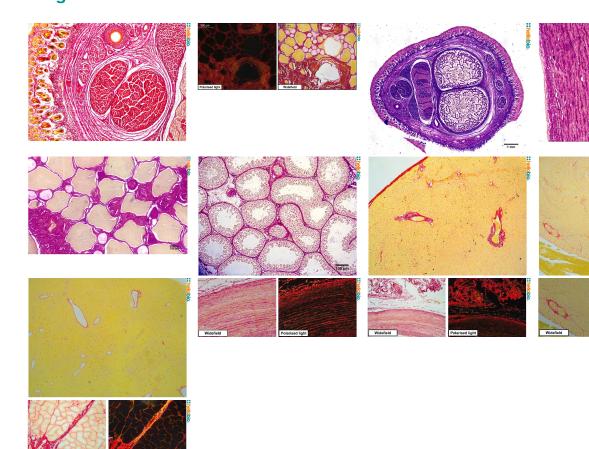

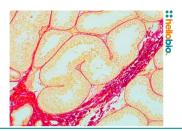

# **Biological Data**

### **Application notes**

#### **Protocol summary:**

- Deparaffinize/ dewax sections if required
- · Hydrate in distilled water
- Immerse sections in picro Sirius red solution and stain for 60 min at room temperature
- Rinse slides quickly in 2 changes of acetic acid solution
- Rinse slide in absolute alcohol
- Deyhydrate in 2 changes of absolute alcohol
- Clear and mount slide using resinous mounting medium. A permanent and not aqueous mounting medium should be used.

#### Interpretation of staining / results:

#### Polarized Light Microscopy:

- Type I (Thick fibers) Yellow-Orange Birefringence
- Type III (Thin fibers) Green Birefringence

### **Light Microscopy:**

- Collagen Red
- Muscle Fibers Yellow
- · Cytoplasm Yellow

Control tissue: Lung, Muscle, Kidney or Uterus may be used as a control tissue.

# **Solubility & Handling**

Storage instructions Important Room temperature

This product is for RESEARCH USE ONLY and is not intended for the rapeutic or diagnostic use. Not

for human or veterinary use

### References

# Method for Picrosirius Red-Polarization Detection of Collagen Fibers in Tissue Sections

Rittie (2017) Methods Mol Biol 1627

**PubMedID** 28836216

# Picrosirius Red Staining: A Useful Tool to Appraise Collagen Networks in Normal and Pathological Tissues

Lattouf et al (2014) J Histochem Cytochem 62(10) **PubMedID**25023614

## Picrosirius Staining Plus Polarization Microscopy, a Specific Method for Collagen Detection in Tissue Sections

Junqueria et al (1979) Histochem J 11(4) **PubMedID**91593